




Aliya® EX + INUMI® Flex for soft tissue ablation

Aliya PEF with the INUMI Flex Endoscopic Needle is designed to integrate into existing endoscopy workflows, enabling biopsy, pathology confirmation, and ablation in the same anesthetic setting.

Delivered through a 21G endoscopic needle, the Aliya PEF system supports predictable ablation near sensitive anatomy while integrating into multidisciplinary care pathways.

Safety & Predictability of Endoscopic Approach


- Compatible with endoscopic tools and workflows
- Minimally invasive endoscopic intervention
- Favorable safety profile, including near sensitive structures
- Preserves future treatment options
- Predictable non-thermal ablation
- Adjunctive component of multimodal care
- Potential for immune activation

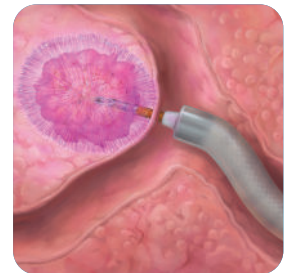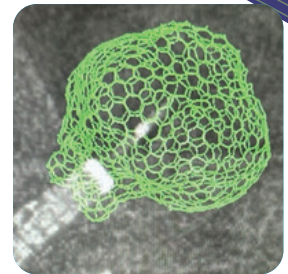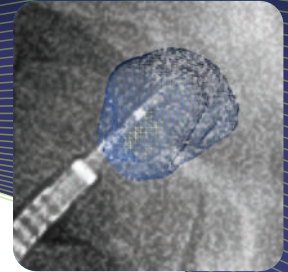


ALIYA EX ABLATION ZONE DIMENSIONS

Approximate PEF volume for one energy delivery based on *in vivo* porcine data.
Representative images, not to scale.


G10 Energy Setting

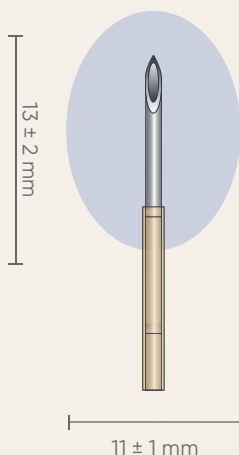

G20 Energy Setting

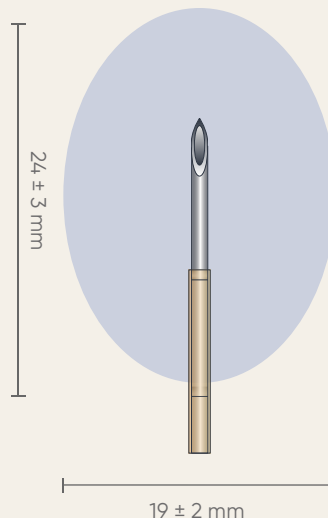

INUMI® FLEX ENDOSCOPIC NEEDLE

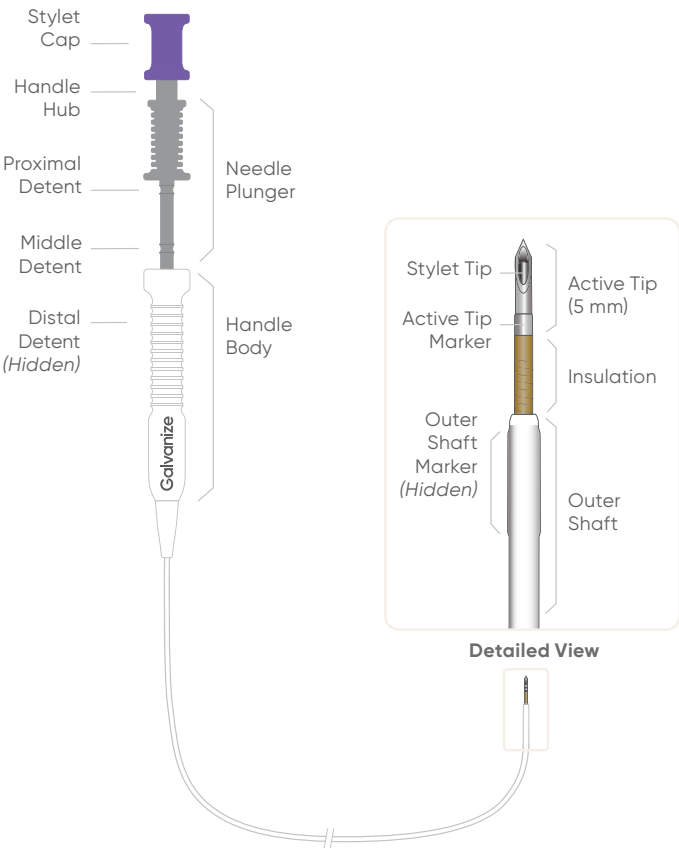

INUMI FLEX NEEDLE

Details	Spec
Needle Gauge	21
Outer Shaft Max Diameter	1.8 mm
Working Channel Min Diameter	2.0 mm
Working Length	160 cm
Active Tip Length	5 mm
Max Activations	8

ALIYA EX ENERGY SETTINGS

Energy Setting	Ablation Dimensions
G10	13 ± 2 mm Length / 11 ± 1 mm Diameter
G20	24 ± 3 mm Length / 19 ± 2 mm Diameter

Approximate PEF volume for one energy delivery based on *in vivo* porcine data.

ORDERING INFORMATION

Part #	Description	UOM
ALIYA EX Starter Pack	Aliya EX Generator	1 Each
	Cardiac Monitor, Synchronization, IVY Model 7600	1 Each
	Galvanize Cart	1 Each
GTI-00054	INUMI Flex Needle	1 Box (5)
GTI-00078	INUMI Flex PEF Ablation Pack: • INUMI Flex Needle, 21G • Aliya Electrode	1 Box (5)

The Aliya® PEF System (K212871, K233884, and K253826) is 510(k) cleared in the United States for the surgical ablation of soft tissue. Galvanize Therapeutics® does not promote the off-label use of its products and nothing herein is intended to promote an off-label use of the Aliya PEF System. The Aliya PEF System is a tool for the surgical ablation of soft tissue it is not intended to treat, cure, prevent, nor mitigate any specific disease or condition.

HANDLE POSITION

